
**VES**electric

# Dental Water Jet / Munddusche Ирригатор полости рта Irrigateur buccal

Model / Modell  
Модель / Modèle: VIP-006-W

EN Instruction manual  
DE Bedienungsanleitung  
RU Инструкция по эксплуатации  
FR Mode d'emploi
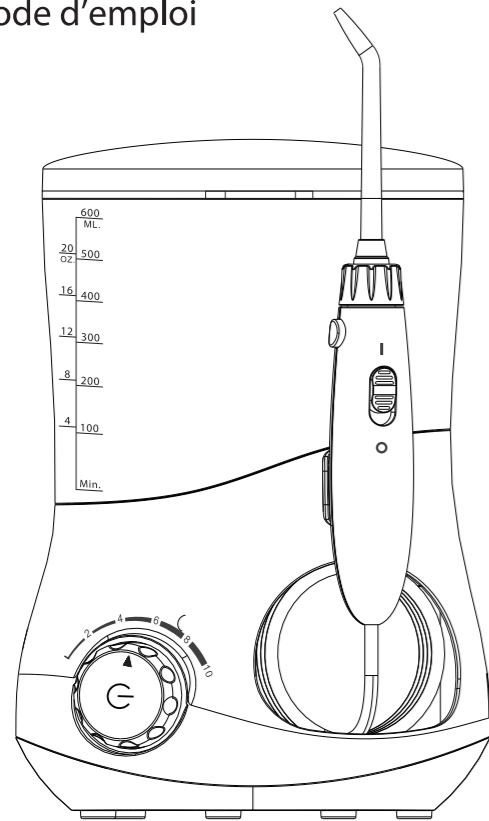
## ENGLISH

### IMPORTANT SAFEGUARDS


Please read the safety instructions before using this product. The safety instructions provided here help you to use this product correctly, and also to prevent you from injuring yourself or others. When using electrical products, especially when children are present, the basic safety precautions should always be followed.

### WARNINGS


Before the first connecting of the appliance to the network check whether the rating indicated on the technical label on the appliance corresponds to power supply in your network.

- For household use only, not for outside use.
- Children should not be allowed to use this device. Please keep it out of reach of children.
- Do not use if the power cord or plug was damaged.
- Always unplug the device after using.
- Do not touch plug with wet hand.
- Do not place the device in water or other liquids.
- Do not place the device where it may fall or be pulled into a tub or sink.
- Do not use the device while bathing.
- Do not touch the product which has fallen into water, Unplug immediately.
- The device can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under supervision.
- Children should not play with the appliance.
- The appliance is only to be used with the power supply unit provided with the appliance.
- Keep the device away from heated surface.
- Do not use other tips except recommend parts.
- Please do not use the device if you lost the nozzle.
- Keep nozzle away from small children. The nozzle is small enough to be swallowed and can be become lodged in the throat.
- Do not direct water under the tongue, into the ear, nose or other delicate areas.
- Do not use for any other purpose except oral cleaning.
- Do not use the appliance: if there is a wound in the mouth; if you wear any items like jewelry on the lips, teeth, mouth remove them before the use.
- People who experience tooth or gum pain may not be able to use this device and should consult their dentist.
- Fill water tank with water or solutions (rinsing or treatment) which are advised by your dentist. The liquid should not be viscous or precipitate otherwise the appliance will be clogged. Do not fill the tank with water which temperature exceeds 40°C. It can cause burns of mucous.
- Do not press the nozzle too close against your teeth or gums, doing so may injure your gums.
- Fill the water tank with filtered water. Tap water contains an admixture of salts, lime, sand, etc., which will quickly clog up the appliance and the nozzle. Appliance, with the dirt from dirty water or mouthwash automatically voids warranty.
- Do not use abrasive, chemical, or alcohol-containing cleaners, polishes to clean the appliance.
- Never attempt to repair the appliance yourself! Warranty will be cancelled automatically if the opening of the appliance is detected. Appliance repair should be performed in authorized service center only!

### PARTS IDENTIFICATION
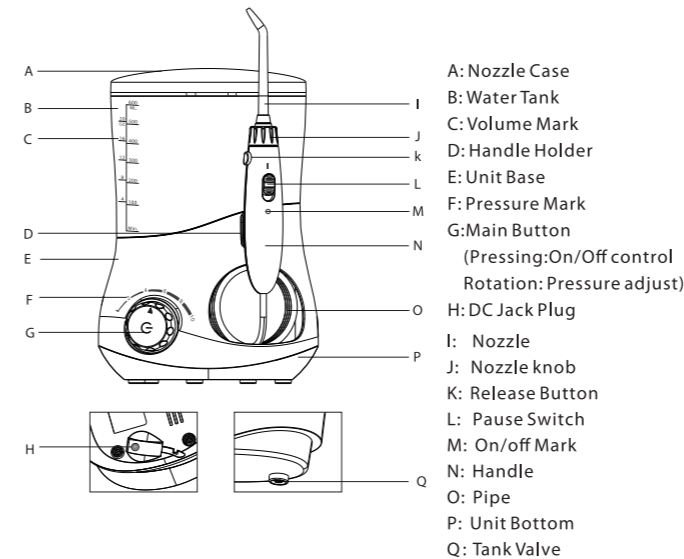
### Accessories
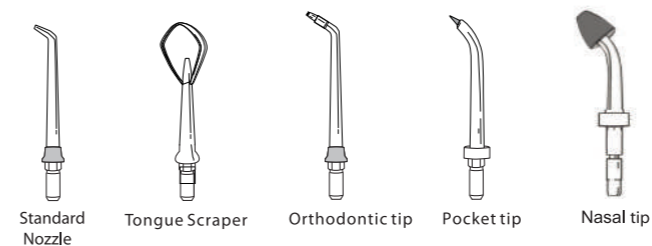
### HOW TO USE

Place the device on a flat surface without any sharp objects and make sure the cord is dry and has no damage. The device is powered by a special adaptor with DC jack connector. (Input:100-240V~,50/60Hz; Output:12V, 15 00mA)

1. Insert the plug into the DC jack on the bottom of base and plug the adaptor cord into an electrical outlet.(Fig.1)
2. Take the Nozzle Case off.
3. Take the Water Tank off to fill with proper liquid, no warmer than 40 °C, do not pour in salt, tooth-paste, or chemicals, then replace the tank onto the top of main unit and press down firmly.(Fig.2-1&Fig2-2)

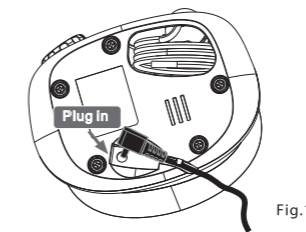
Fig.1

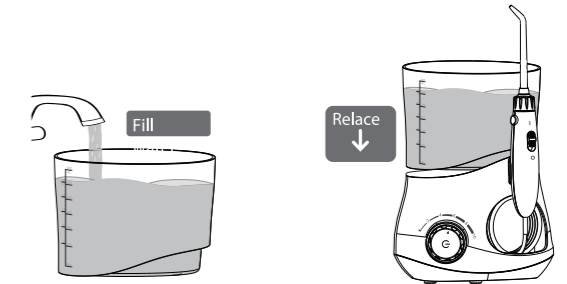
Fig.2-1

Fig.2-2

### Inserting & Removing Nozzle

**Inserting Nozzle**  
Insert tip into the center of the handle, audibly indicating when the nozzle attached to the handle in place. And also try to pull the tip out gently to confirm the tip is locked firmly.(Fig.3)

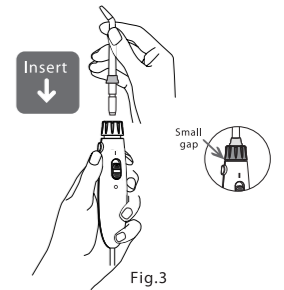
Fig.3

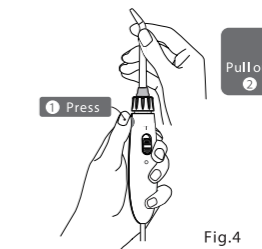
Fig.4

**Removing Nozzle**  
Press the Release Button on the top end of the handle, then pull the nozzle out lightly by hand. (Fig.4)

### How to use the water flosser

**Note:**  
Do not power on the device before placing the nozzle in your mouth. First use, we suggest you choose the lower pressure, gradually increase the water pressure over time to the setting you need.

1. Hold the handle and point the nozzle to sink, then press down the Main Button to power the device on. (Fig.6)
2. Rotate the Main Button to adjust water pressure, after setting, push the Pause Switch to "I" position to start irrigation. To pause the water flow push the Pause Switch to "0" position. (Fig.7-1&Fig.7-2)

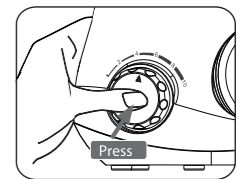
Fig.6

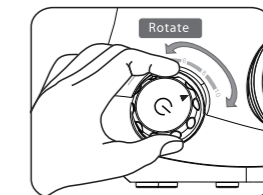
Fig.7-1

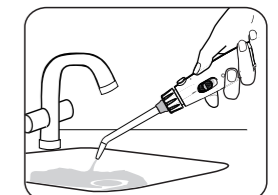
Fig.7-2

3. Adjust the direction of nozzle by rotating the nozzle knob;

Fig.8

4. Bend over to the sink and place the tip in your mouth, and then push the Pause Switch to "I" position to enjoy irrigating. (Slightly close your mouth to avoid the water splashing.) (Fig.8)

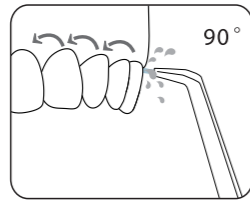
Fig.9

5. Position the water jet flow at 90-degree angle to the teeth and gums. (Fig.9)

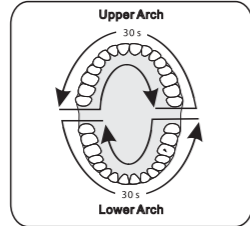
Fig.10

6. Move the water jet slowly along the teeth. The water flow will short pause at **30 seconds**, which reminds you to clean the other arch of your mouth and another short pause at **1-minute** to remind you the recommended usage time is up. The unit will continue to work till **2-minute** cycle is complete or you turn it off. (Fig.10) **Continuous work can last 2 cycles (4 minutes) after which the appliance must cool down for 1 hour.**

#### How to end use

**You can end irrigating through the below methods:**

- 1: Auto shut-off turns off after **2-minute** cycle is complete;
- 2: Turn it off by pressing the Main Button;
- 3: Unplug the adaptor directly when in any emergency;

#### After use

1. Take off the water tank and dispose of water in it.
2. Push the Pause Switch to "I" position, and power on to drain the water inside the device, then push the Pause Switch to "0" position.
3. Turn off the device and wipe it with a soft dry cloth. (Fig.11)
4. Press down the Release Button, and remove the nozzle, then place it into the nozzle case. (Fig.12)
5. Unplug the adaptor.

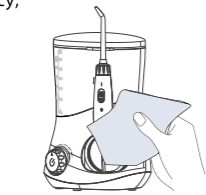
Fig.11

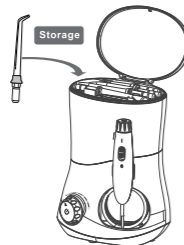
Fig.12

## FUNCTIONAL TIPS

### Tongue Scraper

When you want to clean your tongue, please gently place the tongue scraper onto the center midline of your tongue surface about halfway back with light water pressure setting first, and gradually increase the pressure as needed. (Fig.1 3)

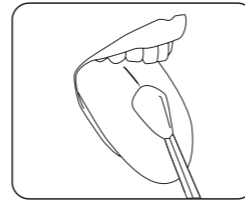
Fig.1 3

### Orthodontic Tip

The orthodontic tip is specially designed for people wearing braces: it will clean the hard to reach areas around and in braces.

- 1: Aim water stream at 90 degree to the tooth at the gum-line, start on the last molar (back tooth);
- 2: Slowly slide the tip along with the gum-line, and pausing briefly to use the brush bristles to clean the area between teeth and all around orthodontic braces.

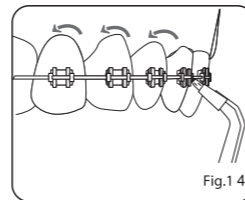
Fig.1 4

- 3: Repeat until you have finished cleaning the inside and outside of all your teeth. (Fig.1 4)

### Tooth Pocket Tip

The pocket tip is specifically designed for the cleaning of periodontal pockets. To use the pocket tip, it is a MUST to set the device to the lowest pressure setting.

1. Place the soft tip against a tooth at 45-degree angle and gently place the tip under the gum-line, then insert into the pocket. (Fig. 15)

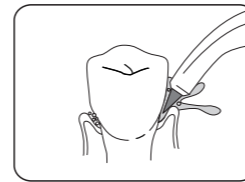
Fig.1 5

2. Turn the unit on and continue tracing along the gum-line until all periodontal pockets are treated. (Fig.1 6)

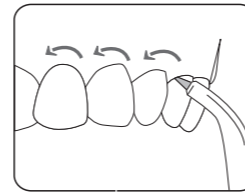
Fig.1 6

### Note:

This tip must be used only with the lowest water pressure to extend life of soft tip as it is easily broken.

### Nasal tip

Nasal tip is used for nose rinsing and ENT diseases treatment. Rinsing nose with salty water or water with medications, unblocks the nose, releases breathe, removes mucus from nasal sinuses. Use it in accordance with below pictures.

## MAINTENANCE

Use only water or neutral cleaning agent to clean the device. Never use caustic or abrasive agents (e.g., vinegar or lime scale removers), since they may damage the device. Do not wash using hot water of over 40 °C. Keep the device away from high temperature area and avoid direct sunlight.

**If you used a mouthwash**, it is necessary to drain the remnants and clean the appliance. To do this, fill the tank with clean filtered water and run water through the used nozzle, including the unit. One hour later, it is necessary to drain clean filtered water through the appliance and the used nozzle.

**Main unit** 1. Remove the water tank and drain the remains water inside the device. 2. Clean by using a soft cloth with liquid soap and dry with a soft dry cloth. 3. Do not submerge the main unit in water to clean.

**Handle** 1. Wash with water and dry up with a soft cloth. 2. Do not bend, pull or twist the water hose.

**Water tank** 1. Wash with water and wipe with a soft dry cloth. 2. Make sure the water tank is inverted and firmly fixed to store. **Nozzle** 1. Wash with water and wipe with a soft dry cloth. 2. Place the nozzle in the storage compartment.

**Charging adaptor** Wipe the metal part of the power plug with a dry cloth once every 6 months.

### Disposal

Dispose the appliance in an environmentally friendly manner. It should not be disposed with the normal household waste. Please dispose it at a recycling centre for used electrical and electronic devices. You can obtain more information from your local authorities.

### Technical data:

Adapter: input: 100-240V, 50/60Hz, 0.8A, output: 12V === 1500 mA  
Max. pressure of water jet: 790 kPa

## TROUBLESHOOTING

Problem	Possible Cause	Solution
Unit does not work	Motor not running	Make sure the adapter is firmly connected to the outlet and the unit.
		Make sure the outlet and unit power is on
Water pressure is very low	Dirt stuck in the nozzle	Clean or replace with a new nozzle
	The nozzle is deformed	Replace with a new nozzle
	Dirt stuck in the tank valve	Clean the valve and fill the tank with pure water
Water does not come out	The tank is empty	Fill the tank up
	Water pressure under the lowest setting	Slide the switch to adjust pressure
	Tank valve not open	Press down the tank, or reassemble the tank valve
Tank leaks	The valve was lost	Replace with a new tank valve
	The valve isn't install in place	Take off the valve and reinstall again
Tear in Pocket tip head	Tip head is worn	Replace with a new nozzle

**NOTE: Manufacturer reserves the right to change parts and color of the appliance. SAVE THIS INSTRUCTION**

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