PORTABLE IRRIGATOR MINI FLOSSER

INSTRUCTION (EN)

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The WATERDENT Mini Flosser provides effective cleaning of the space between teeth, plaque removal and gum pocket cleaning. The using of irrigator prevents dental calculus as well as caries and periodontal diseases.

If you have braces, implants and dental structures, the irrigator Mini Flosser will be an essential tool for oral care. Portable irrigator eliminates harmful bacteria that can cause many diseases of the oral cavity.

Irrigator works from the built-in rechargeable battery. Perfect for both home use and travel.

Please read all instructions carefully before using this device.

- 3. Safety directions
- 5. Indications and contraindications
- 6. Product Descriptions
- 7. Product Accessories
- 8. Before you begin
- 8. Battery charging
- 9. Device usage
- 12. Cleaning guidelines
- 13. Battery maintenance and usage
- 13. Problem solving
- 16. Technical specification
- 17. Symbols explanation
- 17. Warranty statement

Precautionary measures

- Do not use the charger adapter with wet hands
- Do not put the device into the water or any other liquids
- Periodically check the charger adapter cord for damages.
- Do not use the device if there is cord or plug damage
- Keep the charger away from heating devices
- Keep the device away from children
- If your children are under 14, please watch over them during the usage of this device. Give careful instruction of the proper use of the device to your children. It is not recommended using the product for children without a dentist's recommendation
- Persons who cannot use the device independently, as well as persons with a weakened sensitivity of the oral cavity should not use this device
- Use the device only according the given instructions or your dentist recommendations
- Do not put any foreign objects into the recess of the device
- Do not use the device near spray aerosols
- Clean the charger adapter regularly to prevent dust
- Always check whether the charger is connected to a power source with an appropriate rated voltage. Failure to follow this rule may result in electric shock or fire
- Do not disassemble or repair the device by yourself. It may cause a fire, electric shock or injury
- Do not keep the device at places where it may fall or break
- Do not use the device if the jet tips are lost
- Use the device only for cleaning the oral cavity. Water pressure turned to the eyes or ears can cause serious injury
- In case of the oral cavity treatment, as well as severe forms of periodontal disease, please consult your dentist before using this device

- While using do not press the irrigator tip too hard to the teeth or gums, this may cause damage
- When turning off the charger adapter from a power outlet, always hold the power plug (not the power cord)
- Maintain hygienic cleanliness of the device and use only fresh and lightly heated tap water to clean the oral cavity or WATER-DENT irrigator liquids. Some types of liquids can damage the device
- After using the device, be sure to pour out the remaining water from the water tank and turn on the device for 2-3 seconds to complete drying. Failure to follow sanitary-hygienic standards of using the device will result appearance of mildew and bacteria, which can do harm to your health

ATTENTION! As liquids for oral cavity irrigation is strictly **prohibited** to use:

- Substances and liquids containing particulate matters (suspensions, extracts etc.). Such particulate matters can lead to blocking of suction pipe and water pump
- Oily liquids. The oil is an organic solvent, its use may lead to a corrosive effect on the suction pipe and elements of device
- Ordinary and non-matching liquids for mouth washing
- Saline solutions, toothpaste or chemicals
- Mouthwashes that contain iodine

As liquids for the irrigation of the oral cavity is allowed to use:

- Pure warm water. Do not fill the water over 40 ° C into the water tank. In case the quality of tap water in your area is poor, we recommend using filtered water
- WATERDENT irrigator liquids (always mixed with water in a certain proportion. The recommended proportion see the instructions of the mouthwash's manufacturer)

ATTENTION! If you use special liquids for the irrigator, ensure to clean the device with clean warm water after using such liquids. It is necessary to: fill the water tank fully, carry out an additional irrigation procedure, dry out the device with the help of the jet tip. Pay attention to the lifespan of approved for usage liquids. Do not use liquids that have expired.

Indications and contraindications


Please consult your dentist before using the irrigator WATERDENT Mini Flosser.

Irrigator WATERDENT Mini Flosser is to be used for:

- Total care and hygiene of the oral cavity and interdental spaces;
- Care of the oral cavity if you have braces, implants, crowns;
- Gums cleansing and massaging.

Contraindications

- During the first procedures of using the irrigator, insignificant bleeding gum may appear. If you have healthy gums, bleeding gums will stop within 1-2 weeks
- If bleeding continues after 2 weeks of use, stop using the irrigator and consult your dentist
- It is not recommended to use the product for children without a dentist's recommendation, as well as for people with mouth ulcers or open wounds of the oral cavity

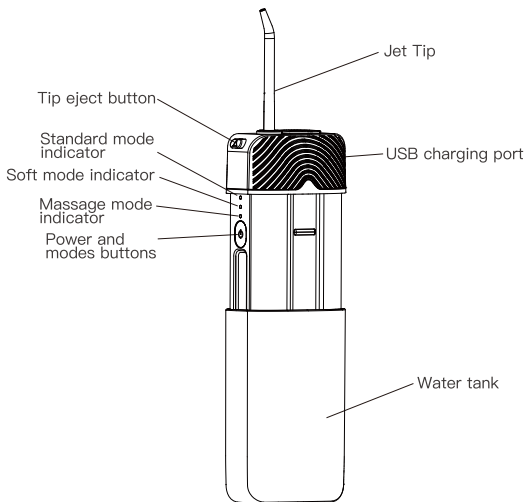
Tip × 2

Charging cable × 1

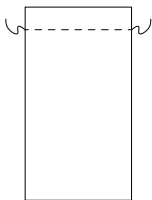Waterproof bag × 1


WATERDENT Mini Flosser includes:

- Irrigator device – 1 pcs.
- Water tank 150 ml. – 1 pcs.
- Jet tip – 2 pcs
- USB-cable– 1 pcs.
- Packaging box – 1 pcs
- Waterproof bag – 1 pcs

What is available:

- On/off switch button
- Modes display light
- Normal mode display light
- Soft mode display light
- Pulse mode display light

Before you begin

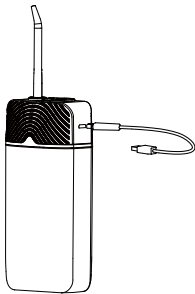Extremely important! Before the first use, please fully charge the battery, insert the USB plug into the charging port, and connect the device to the power source. During charging, the device cannot be used.

Mode Indicator

- Mode indicator. Pressing the on/off switch button will make mode getting change. Chosen mode's indicator will light up
- Normal mode for effective oral cavity cleaning, removal of food pieces stuck between the teeth
- Soft mode for gentle cleaning for sensitive teeth and gums
- Pulse mode for gums cleaning and massaging
- Running in any mode more than 120 seconds, the device **shut down automatically**
- **At any mode press the button more than 1 second, the device will be off**

NB! If the device is used for the first time, to turn it on please press the on / off button, check each of the modes with the filled water tank.

Battery charging


- Connect the USB – cable with charging adapter to the power source

- Connect the device with USB cable
- Mode indicator blinks while the irrigator is charging, when the device is fully charged, indicator will light continuously
- When charging is complete, unplug the charger
- When the device is charging, 3 indicator lights blink. When 3 lights are on, the device is fully charged.

NB! If the indicator lights blink during and after the procedure, it means that the irrigator battery needs to be recharged. It takes 4 hours to fully charge the battery. The battery lasts up to a month of use.

Device usage

Put the jet tip into the correct recess of the main corpus until it clicks. To remove the jet tip, press the tip's eject button and pull up the jet tip

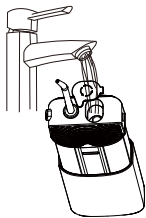Open the lid of the water tank (as shown in the picture) and fill the water or WATERDENT Mouthwash

NB! Please ensure that the device is turned off when filling the water into the water tank

To change or remove a jet tips,
please press the tip's eject button
and remove it.
Fill the water tank with one of the
following ways:

- Open the lid of the water tank. Holding the device in a vertical position, fill the water tank with water of a temperature not exceeding 40 ° C or WATERDENT mouthwash. Close the lid and ensure that it is closed tightly



- Remove the water tank by pulling it. Fill the water tank with water of a temperature not exceeding 40 ° or WATERDENT mouthwash. Put the water tank back into the device

NB! Do not use a water-filled irrigator without a jet tip. This can lead to breakage!

- Before turning on the device, place the irrigator jet tip in the mouth. Bend down over the bathroom sink and keep the device strictly upright, with the jet tip on the teeth and gums. Do not press the on / off button until you insert the tip into the oral cavity

- Start the procedure from the back side of the jaw, stopping at the interdental spaces for better disinfection. Do not stop until all the areas around and between the teeth are cleaned. If you turn the device upside down, water will not leak-in from the water tank. To prevent splashing, during irrigation close the mouth a little

- Direct the water jet as perpendicularly as possible into the interdental space or periodontal pocket. Slowly move along the dentition. Particular attention should be paid to the interdental space and the area around the braces. For the best results start with the molar teeth (posterior teeth), moving gradually to the front teeth. Move along the gums, stopping at the interdental spaces. Do not stop until all the areas around and between the teeth are completely cleaned.

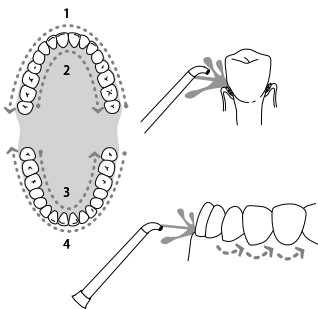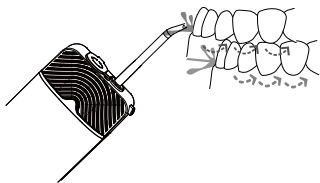- However, remember that mouth irrigation cannot be used instead teeth brushing. Brush your teeth or dentures after a meal for at least three minutes, and then use this device.
- Regardless of the mode, device will automatically turn off within 120 seconds after the using
- After use, press the on / off button to turn off the irrigator. Pull the irrigator out of the mouth
- For ease of use, you can turn the jet tips at an angle of 360

Cleaning guidelines

Main corpus cleaning

1. Turn off the power. Remove the water tank and the jet tips
2. Use only water of a temperature not exceeding 40 ° C or neutral cleaning fluids. It is not recommended to use caustic soda or any abrasive substances to clean the device, as they can damage it
3. Use a piece of cloth to clean the main corpus. Use soap or liquid soap to clean out hard-to-get dirt. Use only neutral cleaning fluids
 - After washing with water, wipe the corpus completely dry with a piece of cloth. If you use any liquids during irrigation, be sure to rinse the device after such use
 - To do this, fill the water tank completely with clean warm water and turn on the device, directing the jet tips down into the bathroom sink, then wipe the appliance dry
4. Rinse the suction pipe with water and wipe with a soft cloth. Do not bend, pull or twist

Water tank

- Rinse water tank with water
- Empty the water tank after cleaning
- If you do not intend to use the device for a week or longer, dry the water tank from the inside and outside
- Use only gentle neutral cleaning substance

Jet tip

- After each use rinse the jet tips under with water
- It is recommended to disinfect jet tips by placing them for 5-10 minutes in one of the following solutions: alcohol solution, 3% hydrogen peroxide solution, 1% chlorhexidine solution
- After disinfection in solution, it is recommended to rinse jet tips with warm water. You can also use non-contact methods of cleaning and disinfecting tips: in an ultrasonic or ultraviolet cleaner

- It is recommended to replace the jet tip with new one every 6 months

NB! When using the irrigator with solutions, be sure to rinse the device after each such use. To do this, fill the water tank completely with clean warm water and turn on the device, at the same time, directing the irrigator tip down to the sink, wipe the device dry.

Battery maintenance and usage

- Before the first use, device must be fully charged (approximately 4 hours)
- If it is necessary to keep the device for a long period of time, please fully charge the device in advance
- To extend battery lifespan, in case the battery is low, please charge the device on time

Problem solving

(EN)

Problem	Reason	Solution
Irrigator does not work	You have just bought a device / it has not been used for a long time (low battery)	Charge the device for 4 hours
Irrigator works only for few minutes even after the charging	Battery lifespan is over	Contact authorized WATERDENT distributor in your country
	Irrigator's battery is low	Charge the device for 4 hours

Water leakage between the jet tip and the cavity for the jet tip	The filler ring on the fixed jet tip cavity is damaged and the tip cannot be installed properly	Contact authorized WATERDENT distributor in your country and install the tip properly
Water leakage of the water tank	Filler ring is weakened	Install properly
Device does not pour out the water	Filtering mesh is plugged up	To clean filtering mesh
Water leakage from the side of the water tank	Inlet port for the air produces water	Common occurrence, the principle of high air pressure
Water pressure extremely low	Jet tip is plugged up	Change jet tip
	Jet tip is deformed	Change jet tip
	Filter is plugged up	Clean out the filter
	The irrigator battery is low	Charge the device for 4 hours

To check the technical condition of the device or in case of a problem that cannot be solved with the help of the abovementioned table, contact an authorized WATERDENT distributor in your country. You can find the address of the authorized WATERDENT distributor on the website <http://waterdent.com/>, as well as on the certificate of guarantee or by telephone +7 495 139 89 13

Technical specification

Charging time	4 hours
Accumulator	Lithium-type battery 1100 mAh
Battery use with full charge	Up to 1 month of usage
Water tank capacity	150 ml
Using time with full filled water tank	Normal mode 40-60 seconds Soft mode 85 seconds Pulse 100 seconds
Water pulse frequency	Up to 1800 imp / min
Max jet pressure	100 PSI

Automatic shutdown timer	After 2 minutes
Condition of use	0°C to 40°C with relative humidity no more than 85%

Symbols explanation

	FDA Certification
	CE Certification
RoHS	RoHS Certification
	FCC Declaration of Conformity
	Electric shock protection class (II class)
	Fragile
	For the duration of shipping/delivery, the product should face upright.
	Keep dry
	Protect from heat
	Special dispose required

Warranty statement

Thank you for purchasing Oral Irrigator WATERDENT Mini Flosser!

Please pay attention to the following aspects:

1. Keep the valid receipt given by the seller as maintenance document

2. This device carries a one-year manufacturer's warranty from the date of purchase (for device malfunction and manufacturing defects)

3. If the product was purchased on the Amazon, no documentation is required. Only the order number is required for us to provide appropriate after-sales service. Service is free of charge (based on the date of the purchase receipt)

4. If the following circumstances arise during the service period, the costs and fees for services are charged depending on the circumstances:

- Damage caused by improper use or maintenance
- Damage caused during self-repair and / or during repair by unauthorized person
- Absence of trade organization's stamp in the warranty card or more than two years have passed from the date of manufacture
- Damage caused by using for the commercial purpose
- Damage caused by force-majeur circumstances (storms, earthquakes, fire)

If irrigator has manufacturing defects, WATERDENT guarantees its free of charge repair or replacement (according to the manufacturer's rules and the type of breakdown) for 12 months from the date of purchase of the device.

The warranty is valid only if the product purchased in Russia or Belarus.

The given warranty coverage does not include:

- Damage caused by improper use
- Normal wear and tear
- Defects that do not significantly affect the use of the device
- Traces of mechanical impact, dimples, flaws, chips on the device
- Damage caused by the impact of the aggressive environment
- This warranty is void if the device was repaired by an unauthorized person or non-original WATERDENT repair parts were used in the repair


If you have certificate of guarantee proving your purchase, contact authorized WATERDENT distributor in your country for warranty maintenance.

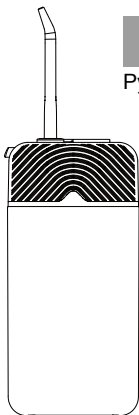All shipping, insurance and any other costs that are relevant to the warranty service will be borne by the buyer. If the case is recognized as a warranty case the buyer is compensated the cost of goods shipping upon the written request of the buyer.

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