**WATER<sup>®</sup>**  
**DENT**


**WATER<sup>®</sup>**  
**DENT**

## SONIC SMART CARE

Sonic electric  
toothbrush

INSTRUCTION (EN)

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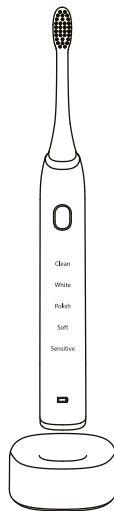
# WATERDENT SONIC SMART CARE

Usage manual

### Thank you for purchasing WATERDENT Sonic electric toothbrush!

Toothbrush WATERDENT Sonic Smart Care with patented DuPont™ bristles and 5 modes of use provides an effective cleaning due to the contact of the bristles with enamel and effect of high-frequency vibrations on dental plaque. Toothbrush prevents tartar and caries formation. Medium bristles indicator shows when brush replacement is required
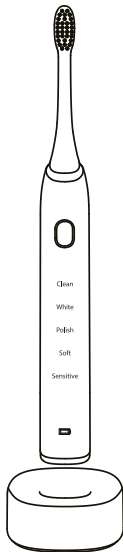
## Contents

Safety directions .....	03
Contraindications.....	04
Product description .....	05
Before you begin .....	06
Battery charging.....	07
Cleaning guidelines.....	09
Battery maintenance and usage .....	10
Technical specification.....	11
Symbols explanation .....	12
Warranty statement .....	12

## Safety directions

- Do not use the charger adapter with wet hands
- Do not put the device into the water or any other liquids
- Do not use while swimming
- Periodically check the charger adapter cord for damages.
- Do not use the device if there is cord or plug damage
- Keep the charger away from heating devices
- Give careful instruction of the proper use of the device to your children. If your children are under 14, please watch over them during the usage of this device
- Persons who cannot use the device independently, as well as persons with a sensitivity of the oral cavity should not use this device
- It is not recommended to use the product for children without a dentist's recommendation, as well as for people with mouth ulcers, open wounds of the oral cavity, in the period after dental surgery or while treatment
- Use the device only according the given instructions or your dentist recommendations
- Do not put any foreign objects into the recess of the device
- Do not use the device near spray aerosols
- Keep the device away from children
- Clean the charger adapter regularly to prevent dust
- Always check whether the charger is connected to a power source with an appropriate rated voltage. Failure to follow this rule may result in electric shock or fire
- Do not disassemble or repair the device by yourself. It may cause a fire, electric shock or injury
- In case of necessary repair contact an authorized WATERDENT distributor in your country
- Do not keep the device at places where it may fall or break
- Do not use the device if the brush heads are lost
- When turning off the charger adaptor from a power outlet, always hold the power plug (not the power cord)

## Safety directions


- This device has only been made for cleaning teeth, gums and tongue. Do not use it for any other purpose
- The toothbrush is a personal care device, suggest to use for use multiple patients in a dental practice or institution
- Stop using a brush head with crushed or bent bristles
- Replace the brush head every 3 months or sooner if signs of bristles fade appear
- Do not use other brush heads instead of the ones recommended by the manufacturer

## Contraindications

Please consult your dentist before using the sonic toothbrush WATERDENT Sonic Smart Care. It is not recommended to use the product:

- For children without a dentist's recommendation, as well as for people with mouth ulcers or open wounds of the oral cavity, in the period after dental surgery or while treatment
- If while using sonic toothbrush bleeding appears, stop using the device and consult your dentist

## Product description
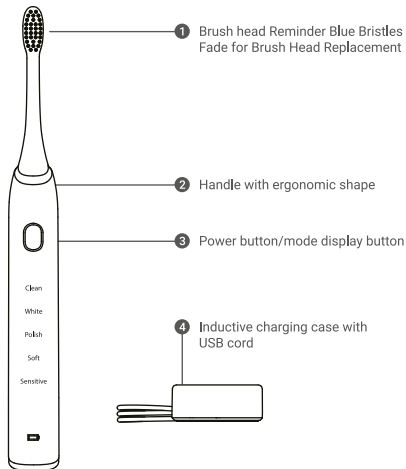
## Before you begin

**EXTREMELY IMPORTANT!** Before the first use, please fully charge the battery, insert the USB plug into the charging port, and connect the device to the power source. During charging, the device cannot be used

## Instruction for use

**01**

Put brush head into toothbrush body and make a certain gap between head and body

**02**

Take some toothpaste on the wet bristle

**03**

Smudge the toothpaste to teeth before use (Notice: launch toothbrush after put into mouth)

**04**

Change the cleaning area based on briefly pause during the brushing cycle (Notice: no need to left and right or upper and down to move after launched)

## Mode Indicator

Pressing the on/off switch button will make mode getting change. Chosen mode's indicator will light up.

- **Clean** - everyday teeth cleaning
- **White** - removes soft plaque
- **Sensitive** – sensitive teeth care
- **Soft** - gums massage
- **Polish** – soft enamel polishing


**NB!** If the device is used for the first time, to turn it on please press the on / off button, check each of the modes

## Battery charging

**NB!** If the indicator lights blink in red color during and after the procedure, it means that the device battery needs to be recharged. It takes 8 hours to fully charge the battery. The battery lasts up to a month of use.

## Battery charging

- Charge for above 12hrs before first use. Please understand that wireless charging is comparably slow for the battery charging.




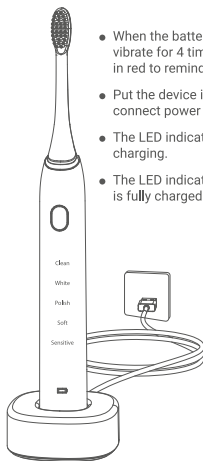
- When the battery power is low, the device will vibrate for 4 times and the LED indicator will flash in red to remind the user to charge the battery.
- Put the device into charging base and confirm connect power source correctly.
- The LED indicator flash in white during battery charging.
- The LED indicator light in white color when battery is fully charged.

Attention:

1. Don't charge the Toothbrush by computer.
2. Don't charge the Toothbrush on Metal desktop.
3. Guarantee the charging time is around 8 hours every time.







## Cleaning guidelines

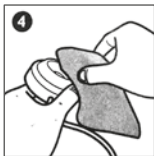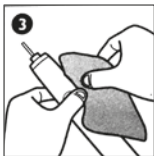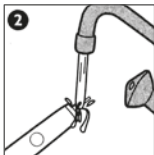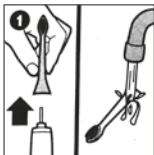Clean the Smart toothbrush periodically. Rinse the handle and the brush head. Wipe down all surfaces with a soft cloth and mild soap, if desired. Do not use sharp objects to clean. Do not clean the brush head, handle, charging base in the dishwasher.

## Battery maintenance and usage

- Before the first use, device must be fully charged (approximately 12 hours)
- If it is necessary to keep the device for a long period of time, please fully charge the device in advance
- To extend battery lifespan, in case the battery is low, please charge the device on time

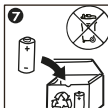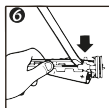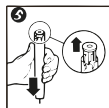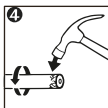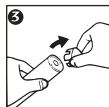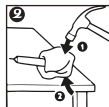Do not throw away the device with the normal household waste at the end of life, but hand it in at an official collection point for recycling. Always remove the battery before you discard and hand in the appliance at an official collection point. This process is not reversible.

To remove the recharge battery, you need a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure above. Be sure to protect your eyes, hands, fingers and the surface on which you work.

## Technical specification

Waterproof	IPX7
Function Operation	Button Press
Outlook	Main body material: PC + ABS + UV Painting
Product Size	Main device : 27.8 * 27.8 * 248.6mm Charging base: 65.6 * 65.6 * 26.3mm
Product Weight	Main body about 115g Charging Base:100g
Charging Type	Wireless Charger
Charging Time	About 8 hours
Battery	Lithium rechargeable Battery
Capacity	800mAh
Battery Life	About 1 month (if use 4 minutes everyday)

Brush Modes	<b>5 different modes, each mode will be lighted up when relevant mode function is on</b>
	Clean mode: 34800 brush strokes per minute White mode: 33000-40000 brush strokes per minute Polish Mode: 35000 brush strokes per minute Soft mode: 33000 brush strokes per minute Sensitive mode: 28200 brush strokes per minute

Package Contents	2 Brush head 1 Handle with ergonomic shape 1 Wireless Charger 1 User Manual 1 Gift Box
------------------	--

## Symbols explanation

	FDA Certification
	CE Certification
<b>RoHS</b>	RoHS Certification
	Eurasian Conformity
	Fragile
	For the duration of shipping/delivery, the product should face upright.
	Keep dry
	Protect from heat
	Special dispose required

## Warranty statement

Thank you for purchasing Sonic electric toothbrush WATERDENT Sonic Smart Care! Please pay attention to the following aspects:

1. Keep the valid receipt given by the seller as maintenance document
2. This device carries a one-year manufacturer's warranty from the date of purchase (for device malfunction and manufacturing defects)

## Warranty statement

3. If the product was purchased on the Amazon, no documentation is required. Only the order number is required for us to provide appropriate after-sales service. Service is free of charge (based on the date of the purchase receipt)

4. If the following circumstances arise during the service period, the costs and fees for services are charged depending on the circumstances:

- Damage caused by improper use or maintenance
- Damage caused during self-repair and / or during repair by unauthorized person
- Absence of trade organization's stamp in the warranty card or more than two years have passed from the date of manufacture
- Damage caused by using for the commercial purpose
- Damage caused by force-majeure circumstances (storms, earthquakes, fire)

If device has manufacturing defects, WATERDENT guarantees its free of charge repair or replacement (according to the manufacturer's rules and the type of breakdown) for 12 months from the date of purchase of the device. The given warranty coverage does not include

- Damage caused by improper use
- Normal wear and tear
- Defects that do not significantly affect the use of the device
- Traces of mechanical impact, dimples, flaws, chips on the device
- Damage caused by the impact of the aggressive environment
- This warranty is void if the device was repaired by an unauthorized person or non-original WATERDENT repair parts were used in the repair

## Warranty statement

If you have certificate of guarantee proving your purchase, contact authorized WATERDENT distributor in your country for warranty maintenance.

All shipping, insurance and any other costs that are relevant to the warranty service will be borne by the buyer. If the case is recognized as a warranty case the buyer is compensated the cost of goods shipping upon the written request of the buyer.

The production date is included in the serial number of the device.

WATERDENT distributor in Russia:  
LLC «Global Dent», Russia, 111524, Moscow, Elektrodnyaya street,  
9B, floor 1, room. 113.  
Tel: +7(495) 139-89-13

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