**EN** Toothbrush **AENO DB7/DB8** is designed for daily cleaning of teeth and gums.

### Specifications

Battery: Non-removable 3.7 V lithium-ion (DC), 800 mAh. Full charge time: 4 h. Power: 0.5 W. Noise level (at 0.3 m): ≤60 dB. Body protection degree: IPX7. Stiffness of the bristles: soft. Number of operation modes: 3. Dimensions of the brush with head: Ø 21 mm × 208 mm.

### Scope of Delivery (figure a, d, e)

Toothbrush **AENO DB7/DB8** (white) (handle (1), removable head (2) (DB7 – 1 pc., DB8 – 3 pcs.)), power cord micro-USB / USB Type-A (0.3 m), quick start guide, warranty card, (DB7 – 2 stickers, DB8 – metal tool (1 pc.), mirror (1 pc.)).

### Limitations and Warnings

This appliance is not intended for use by children or persons with reduced physical, sensory or mental capabilities if they have inadequate experience or knowledge of how to use the appliance, or if they are not under the supervision of a person responsible for their safety. If you have diseases of the teeth, gums and mouth, you should first consult a dentist about the possibility of using the device. Replace the brush head at least every three months, depending on wear. Do not perform maintenance on the device and its components until the device is turned off and disconnected from the power supply. Do not wash the unit or any part of it in the dishwasher. When washing, do not immerse the brush handle and the head completely in water (figure c), do not use aggressive cleaning agents. It is acceptable to occasionally treat the bristles of the head with mouthwash. Dirt on the handle of the device should be removed with a damp cloth. Do not disassemble or attempt to repair the unit yourself.

### Preparing the device for operation

Remove the cover (3). Connect the device to the mains using the USB cable and the AC adapter. When the LED indicator (4) lights green, the battery is fully charged, and the device is ready for use. Charge the device for at least 24 hours before using it for the first time. If you brush your teeth twice a day for 2 minutes, the device will work without a charge for up to 3 months.

**ATTENTION!** The AC adapter is not included in the package. You can use any adapter with a 5.0 V / 0.3 A (DC) power setting.

### Operating the device

Install the replacement head (2) on the handle (1) as shown in figure b. Press the power button (8). The indicator (4) will light up, as well as one of the indicators that show the operating mode of the device. To turn the device off, press the button (8) again.


Every 30 seconds, the vibration alarm goes off, reminding you to switch the operating mode. To change the mode, switch the device off and on again with button (8). Then quickly press the button (8) several times in a row to select the desired mode. Each time you press it, the mode changes to the next.

### Operating modes of the device (figure a)

- "Soft" (5) – delicate cleaning of sensitive areas
- "Clean" (6) – daily standard cleaning
- "White" (7) – teeth whitening

### Troubleshooting

- The replaceable head vibrates poorly. Possible reason: the head fits too tightly on the handle of the device. Solution: make sure that there is a small gap at the mating point (9), about 0.5 mm.
- Short operating time after a full charge. Possible reason: the device has not been used for more than 6 months, so the battery does not hold a charge very well. Solution: contact the service center and replace the device battery.

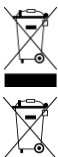These symbols mean that you must follow the Waste Electrical and Electronic Equipment (WEEE) and Waste Battery and Battery Regulations when disposing of the device, its batteries and accumulators, and its electrical and electronic accessories. According to the rules, this equipment must be disposed of separately at the end of its service life.

Do not dispose of the device, its batteries and accumulators, or its electrical and electronic accessories together with unsorted municipal waste, as this will harm the environment. To dispose of this equipment, it must be returned to the point of sale or taken to a local recycling facility. You should contact your local household waste disposal service for details.

**The warranty period and service life is 2 years from the date of purchase of the product.**

Manufacturer information: ASBISc Enterprises PLC, Agios Athanasios, Lapetou 1, 4101 Limassol, Cyprus. Made in China. Up-to-date information and a detailed description of the device, as well as connection instructions, certificates, information on companies that accept quality claims and

warranties, are available for download at [aeno.com/documents](http://aeno.com/documents). All trademarks and brand names are the property of their respective owners.

ASBISc reserves the right to modify the device and to make edits and changes to this document without prior notice to users and shall not be liable for possible negative consequences resulting from the use of an outdated version of the document or for any possible technical or typographical errors or omissions or incidental or consequential damages that may arise from the transmission of this document or the use of devices.

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