
Overview

1. Flat Head (Non-removable)
2. Battery Indicator
3. Charging Port
4. Power Button (Selecting gears/Power switch)
5. Gear Indicator (Breathing light)

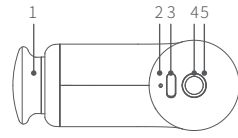Flat Head

Massage flat major muscle groups on the back.

Spherical Head

Relax soft tissues and sensitive muscles.

Cone Head

Percuss deep tissues like meridians and acupoint, etc.

U-shaped Head

Massage neck (not cervical vertebrae) and tendon.

List

- | | |
|-----------------------------------|----|
| ① Main Body (including flat head) | ×1 |
| ② Spherical Head | ×1 |
| ③ Cone Head | ×1 |
| ④ U-shaped Head | ×1 |
| ⑤ USB-C Cable | ×1 |
| ⑥ Manual | ×1 |


Contraindications

1. Cancers.
2. Metal bone anastomosis (do not use the device directly in the fracture area).
3. Osteoporosis.
4. Vein thrombosis and pulsation.
5. Hypertension.
6. Pregnant woman.
7. Varicose veins.
8. Acute inflammation.
9. Do not massage the bony area.
10. Do not massage surgical site.
11. Do not massage the cervical vertebrae and spine.
12. Kidney stones (do not use the device near the kidneys).
13. Do not massage one episode for more than 5 minutes.

Start to Use


1. Massage Head Combination

The device can be used directly with the built-in massage head (flat head) on it, or combine with another one depending on your needs. See complete combination as below.

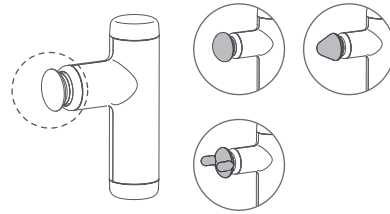2. Power ON/OFF

Power on: Long press the power button.
Power off: Long press power button or single tap it to switch modes until the product power off.

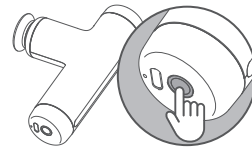3. Selecting Gears

Single tap power button to select gears after powering on.

| Gears | Rated rotate speed (rpm) | Function |
|-------------|--------------------------|---------------------|
| First gear | 2000 | Muscle arousal |
| Second gear | 2400 | Relax |
| Third gear | 3000 | Deeply massage |
| Fourth gear | Auto change | Intelligent massage |
| Power off | 0 | |

| Gears | Gear indicator (breathing light) |
|-------------|----------------------------------|
| First gear | Blue breathing light |
| Second gear | Green breathing light |
| Third gear | Red breathing light |
| Fourth gear | Yellow breathing light |

4. Intelligent Protection

This product has intelligent protection function. In order to prevent muscle injury caused by massaging for a long time, gear indicator will turn off (battery indicator remain) after working for 10 minutes continuously, and it can be restarted by pressing the function button. If there is no operation within 3 seconds, it will power off automatically.

⚠ Do not press the massage head firmly to block the motor rotation while using, otherwise it will reduce the service life of the product!

5. Battery Indicator

| Indicator while using | Description |
|-----------------------------|-------------------------|
| Constant green light | Fully charged |
| Green light flashes slowly | Medium level of battery |
| Green light flashes quickly | Low power |

| Indicator while charging | Description |
|--------------------------|------------------|
| Constant red light | Charging |
| Constant green light | Charge completed |

6. Charging

⚠ Do not use it while charging.

7. Specifications

Product Name: KiCA Massage Device
Product Model: MINI-C
Operating Voltage Range: 7.1V~8.4V
Operating Temperature Range: 0°C~+40°C
Charging Temperature Range: 0°C~+40°C
Battery Capacity: 1300mAh
Charging Time: >2 Hours
Weight: 213g
Size: 88.2×35×128mm
Theoretical Battery Life: About 7 Hours (Test in lab environment under the first gear with flat head)

Safety Instruction

1. It is recommend to massage your body after wearing dry clothes or placing a clean and dry towel on massage part, then gently press and move the device.
2. High-speed and high-pressure massage can only be used to soft tissues without feeling pain or maladjustments. Do not massage head and any skeletal sites with high speed and high pressure.
3. Massage the same part of the body with high speed and high pressure for a long time will.
4. Please keep your fingers, hair or, any other body hair away from the extendable part of the device to prevent pinching and clipping.
5. Do not discard or misuse this device.
6. Do not modify this device by yourself without permission.
7. Do not operate or charge the device when nobody was around.
8. This device is for adults only.
9. People with contraindicated symptoms should use this device with doctor's permission.

Notice

1. Make sure the massage head is not blocked by external force when the device is power on. If it is blocked, it will pause for a while and try to restart.
2. DO NOT contact the product with water or other liquid for the massage device is not waterproof.
3. DO NOT disassemble the product except marked detachable. It need send to after-sales or authorized service center to fix it if you accidentally disassemble it and cause abnormal work. The relevant costs are borne by users.
4. Prolonged continuous operation may cause the motor surface temperature to rise, please operate carefully.
5. DO NOT drop or strike the product. If the product is abnormal, contact After-sales support.

Storage and Maintenance

1. Keep the product out of the reach of children and pets.
2. DO NOT leave the product near heat sources such as a furnace or heater. DO NOT leave the product inside of a vehicle on hot days.
3. Please storage the product in a dry environment.
4. DO NOT overcharge or overuse the battery, otherwise, it will cause damage to the battery core.
5. Never use the product when the temperature is too high or too low.

Precautions when discarding the product

The lithium ion battery used by the product is a recyclable and renewable valuable resource. When discarding it, please have it handled by a professional agency that can recycle chargeable batteries. Before discarding the product, power it off and take the lithium ion battery out from it following the steps below. The battery is to be disposed of safely.

1. Disassemble the top housing with help of a professional tool.
2. Pull out the batter terminals from the PCBA.
3. Use a professional cutter to cut off the bottom housing and then take out the battery.

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Warranty:
(1) Within one year from the date of selling, the product is malfunctioning under normal due to non-artificial reasons.
(2) The malfunction of the product is not caused by artificial reasons such as unauthorized disassemble or conversion or addition.
(3) The buyer is able to provide the certificate of maintenance service: the warranty card, the legitimate receipts, invoices, or screenshot of purchasing.

The followings cases are not included under the warranty:
(1) Unable to provide legitimate receipt and warranty card with buyer's info.
(2) The damage caused by human or irresistible factor.
More details about after-sales policy, please refer to after-sales page on website: <https://www.kica-tech.com/en/support>

* Our company reserves the right of final interpretation on the above-mentioned after sales terms and limitations.

| | | |
|------------------------------|-------------------------------------|------------|
| 合格证 QUALIFIED CERTIFICATE | 检验员 Checker 生产日期 Production Date | QC PASS |
|------------------------------|-------------------------------------|------------|