
# Naipo

## User Manual

Shiatsu 3D Rotating Massager  
with Rechargeable Battery

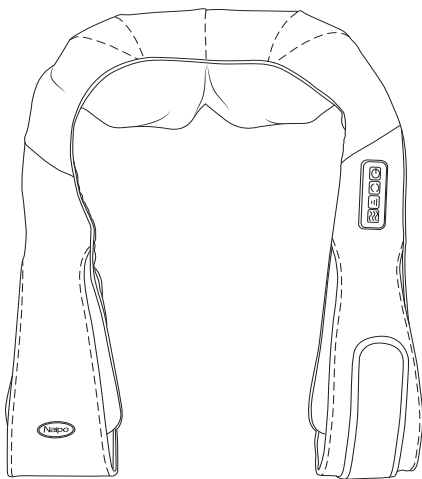Model: MGS-321

# Contents

---

English	01~04
Deutsch	05~08
Français	09~12
Español	13~16
Italiano	17~20

---

**Please make a full charging to it before the first using, because the battery may not have enough power after a long storage.**

## **Important:**

Please keep this manual properly for future reference. Before using this appliance, please read all the instructions carefully.

## **Safety Instructions:**

Please follow basic safety instructions as below to reduce hazards or risks of electrical shock, burns or any other injuries.

- Always unplug the appliance from the electrical outlet immediately after using and before cleaning.
- Do not reach for the appliance which has fallen into water, unplug it immediately.
- Do not use it while bathing or in the shower.
- Do not place or store the appliance where it may fall into a tub or sink.
- Never use pins or other metallic fasteners with this appliance.
- Do not use it in wet or moist condition.
- An appliance should never be left unattended when plugged in. Unplug from outlet when not in use and before putting or taking off parts and attachments.
- Close supervision is necessary when this appliance is used by, on or near children, invalids or disabled persons.
- Use this appliance only for its intended use as described in this manual. Only use attachments provided with the unit.
- Never operate this appliance if it has a damaged cord or plug, or not working properly.
- Keep cord away from heated surfaces.
- Never use it while sleeping or fall asleep while using the massager.
- Never drop or insert any object into any opening.
- Do not operate while aerosol (spray) products are being used or where oxygenize being administered.

- Do not operate under a blanket or pillow. Excessive heating can occur or cause fire, electric shock or injury to persons.
- Do not carry this appliance by power cord or use cord as handle.
- To disconnect, turn all controls to the "OFF" position, then remove plug from outlet.
- This appliance is designed for personal, non-professional use only.
- Do not use outdoors.
- Do not use on an infant or invalid or on a sleeping or unconscious person.
- Do not use it on sensitive skin or somebody who is insensitive to heat.
- Do not pull or twist the appliance on the cable, the power supply unit or the control.

## **Cautions:**

Please read all the instructions carefully before operating.

- This product is not designed for people ( including children ) with limited physical, sensory or mental abilities, unless they are supervised by an adult who responsible for their safety.
- The unit can be used with rechargeable battery only or directly with power adapter.
- If you are experiencing health complication or conditions as stated below, please consult your doctor before use:
  - you are pregnant.
  - you have a cardiac pacemaker, artificial joints or electronic implants.
  - you are suffering one or more of the following illnesses or complaints: circulation disorders, varicose vein, open wounds, bruises, broken skin, inflammation of the veins.
- Do not use this unit to replace medical treatment.
- If you experience any pain or discomfort during the massage, please stop immediately and consult your doctor.
- Never use this product directly on swollen or inflamed areas or skin eruptions.
- Do not use this product before bed. The massage has a stimulating effect and may delay sleep.
- Children or invalids should use this product under the adult supervision.

## Technical Parameters:

Description	Neck & Shoulder Shiatsu Massager with Rechargeable Battery		
Model NO.	MGS-321	Power	25.6W
Input voltage	12.8V DC 2A		
Li-on Battery	11.1V 2000mAh		

## Controller Diagram

Power Button: On / Off

Rotate Direction: Anticlockwise / Clockwise

Massage Speed: Medium / High / Low

Heating Button: On / Off

## Operation Guide

### 1. Using it while charging

- ① Connect the power adapter and the massager correctly, and then plug it into the power supply.
- ② Press the power button to turn on the device, and the default mode is middle rotating speed and heating function. The Lithium battery is not charging while using it.
- ③ Customize the massage functions by selecting rotate direction, speed or heating buttons respectively.
- ④ Press the power button again to turn off the device.

### 2. Using it with rechargeable battery

- ① Long press the power button for 2 seconds to turn on the device, and the default mode is middle rotating speed and heating function.
- ② Customize the massage functions by selecting rotate direction, speed or heating buttons respectively.
- ③ Press the power button again to turn off the device.

\* Full charged Lithium battery can be used continuously for 6 periods. ( 15 mins / per period )


Please make a full charging to it before the first using, because the battery may not have enough power after a long storage. When charging, the blue light turns on, and it will turn off after full charged.

### Massage Areas:


This massager can be applied to use in neck, shoulders, back, waist, abdomen, thighs or legs.


Neck & Shoulder

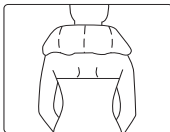
Neck

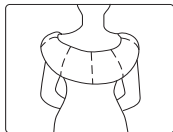Back

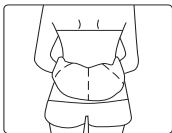Waist

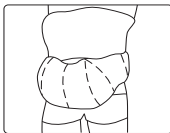Abdomen

Thighs

### Cleaning and maintenance:

Make sure the product is switched off before cleaning. This appliance can be cleaned only by clean soft towel. NEVER use brushes, cleaning agents, petrol, thinners or alcohol. Do not kink the cable.

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Shenzhen Pinlong E-Commerce Co., Ltd.