SCARLETT®

RUS

ENG

KAZ

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The safety measures and instructions contained in this manual do not cover all possible situations that may arise during the operation of the device. When working with the device, the user should be guided by common sense, be careful and attentive.

SECURITY MEASURES

- The manufacturer is not responsible for damages caused by non-compliance with safety regulations and rules of operation of the device.
- This device is designed for domestic use and can be used in apartments, country houses or other similar conditions of non-industrial operation. Industrial or any other inappropriate use of the device will be considered a violation of the conditions of proper operation of the device. In this case, the manufacturer is not responsible for the possible consequences.
- Do not extend the USB cable in doorways or near heat sources. Make sure that the cable does not twist or bend, does not come into contact with sharp objects, corners and edges of furniture.
- Improper handling may result in damage to the device, damage to property and cause harm to the health of the user.
- If the power cord is damaged, it should be replaced with a special cord obtained from the manufacturer or service center. Also the same type of cable can be bought in the shop.
- The manufacturer reserves the right to make minor changes in the design of the device that do not drastically affect its safety, performance and functionality.
- Do not attempt to repair the device yourself. In case of malfunctions, please contact your dealer.
- The date of manufacture is indicated on the device and/or on the packaging, as well as in the accompanying documentation, in the format XX.XXXX, where the first two digits "XX" are the month of manufacture, the next four digits "XXXX". XX" - the year of manufacture.

ATTENTION! THE WATER FLOSSER IS NOT A SUBSTITUTE FOR A TOOTHBRUSH, THEREFORE, IN ORDER TO ENSURE COMPREHENSIVE ORAL CARE OF THE ORAL CAVITY, IT IS RECOMMENDED TO MAKE BRUSHING FIRST. THE DEVICE IS DESIGNED TO CLEAN TEETH AND ORAL CAVITY, USE IT ONLY FOR ITS INTENDED PURPOSE.

REMEMBER: accidental damage to the USB cable can lead to problems that do not comply with the warranty conditions, as well as can cause electric shock.

- It is forbidden to operate the device outdoors: moisture or foreign objects entering the device body may cause serious damage to it. Do not charge the device near water-filled containers (bathtubs, sinks, etc.). Keep the USB cable away from water.

ATTENTION!

- It is not recommended to immerse the device body in water or other liquids for a long time. Before charging, make sure that the device is completely dry.
- Consult a dentist before using the device.
- When wearing dentures, implants, crowns, etc., make sure that they can be cleaned with a water flosser. The manufacturer is not responsible for possible damage.
- Do not direct a stream of water under the tongue, into the ears, nose, eyes and other sensitive areas of the body – this can lead to injury.
- Take out all removable appliances from the mouth before using the device.

- Before cleaning the device, make sure it is unplugged and turned off. Strictly follow the cleaning instructions.
- The device is not intended for use by people (including children) with reduced physical, sensory or mental abilities or with lack of experience or knowledge, unless these people are supervised or instructed in how to use this device by a person responsible for their safety.
- Children should be supervised to prevent playing with the device.
- Self-repair of the device or modification of its design is prohibited. Unprofessionally performed work can lead to damage to the device, injury and property damage.

ATTENTION! IT IS FORBIDDEN TO USE THE DEVICE IN CASE OF ANY MALFUNCTION.

Technical Specifications:

Type of the device: Water Flosser

Model: SC-WF01

Device performance: 3.7V  7.5W

Charging performance: 5V  2W

Mode: 4 modes

Waterproof rate: IPX7

Charge time: ≤ 6 hours

Water pressure range: 552-1103 kPa


Water tank capacity: 350 ml

Overall dimensions (L x W x H): 86 x 81 x 230 mm


Net weight: 307 g

Package contents:

Water flosser	1
Type-C cable	1
Standard nozzle	1
Tongue nozzle	1
Periodontal nozzle	1
Orthodontic nozzle	1
Operation Manual	1
Warranty card	1
Case	1

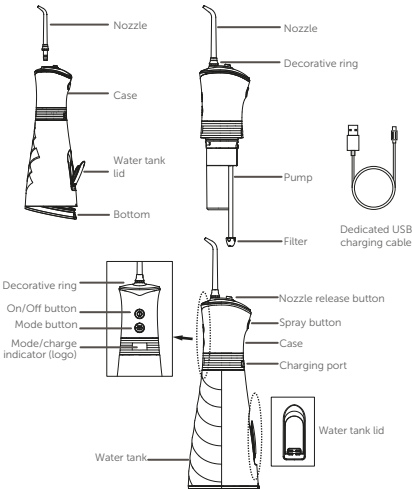Nozzle types:

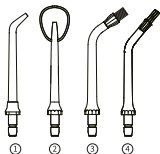
- ① Standard nozzle:
Suitable for teeth cleaning for the majority of people;
- ② Tongue nozzle:
Suitable for cleaning the tongue surface;
- ③ Parodontic nozzle:
Suitable for cleaning dental plaque;
- ④ Orthodontic nozzle:
Suitable for people undergoing orthodontic treatment;

Operation Manual


Charge the device before use as follows:

1. Check and ensure that the device and hands are dry, then connect the USB charging cable to the power adapter and plug the power plug into the power socket. (Figure1)
2. Insert the charging plug of the USB charging cable into the charging port of the teeth punch, and ensure that it is inserted properly. (Figure 2)

Attention! The device will not operate while it is connected to the power supply.

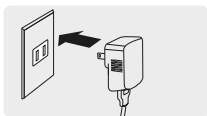(Fig1)

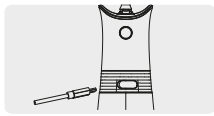(Fig2)

Reminder:

When the red indicator light flashes or the device cannot be turned on properly after use, it indicates that the battery level is low and the device needs to be charged.


·It takes 5-6 hours to fully charge.

·When charging, the red indicator light is on.

·When charging is complete, the green indicator light is on.

·The battery has overcharge protection, so charging will be disabled when the battery is fully charged.

3. When charging is over and the device is unplugged, insert the nozzle into the top hole of the handle (Figure 3-1). If you need to replace the nozzle, turn the product off first and then lightly press the nozzle release button to easily replace it (Figure 3-2).


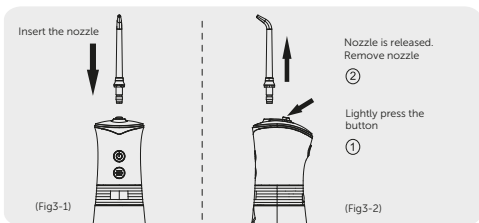
4. Fill the tank with water the following ways:

The first method: Open the tank lid, hold the product horizontally, fill the tank with water and close the lid. (Figure 4)

The second method: Remove the water tank, hold it vertically, and fill it with water. (Figure 5)

ATTENTION! IT IS PROHIBITED TO FILL THE WATER TANK WITH HOT WATER, HERBAL INFUSIONS, SALT, OIL-CONTAINING OR ANTISEPTIC SOLUTIONS, TOOTHPASTE, CHEMICALS, NON-SPECIALISED MOUTH[1]WASHES, AS WELL AS OTHER LIQUIDS NOT INTENDED FOR IRRIGATORS!

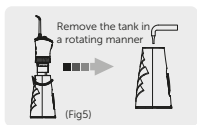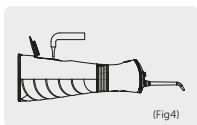5. Turn on the device by pressing On/Off button. In standby mode press the Mode button to choose the desired mode. After 3 seconds of inactivity, standby will be restored.

This water flosser has four modes:

Clean (standard mode): When the mode button is pressed, the logo lights up white (first mode); Standard water pressure is used, suitable for most people.

Soft (soft mode): When mode button is pressed, the logo lights up blue (second mode). Low water pressure is used, suitable for beginners and people with sensitive gums.

Gentle water pressure, suitable for first-time experience and sensitive individuals.

Pulse (pulse mode): When mode button is pressed, the logo flashes blue (third mode). Used for deep cleaning and gum massage.

DIY (manual mode): When mode button is pressed, the logo lights up blue (fourth mode). After the long pressing the mode button the one can adjust the DIY speed, which will be automatically remembered when the button is released.

Pause mode: In operation mode, press and hold the spray button to stop the device from working. Release the spray button to start the device. The pause time is not taken into account when selecting the mode switching time.

The spray button is also designed to spray water briefly. To use, simply press the button once to activate the mode.



Reminder!

For first-time users, use the soft mode for one week first, and then use other modes after adapting.

6. Before use, hold the flosser vertically, open the mouth, align the nozzle with the gaps between the teeth or gums and slightly bend the head so that the water flow out smoothly (Figure 6-1&2).

7. Press the button mode. When using, pay attention to the following points:


- (a) Align the nozzle with the gum line and perpendicular to the gum line;
- (b) Slowly move along the teeth and clean them in order;
- (c) Thoroughly clean teeth, braces, crowns, and dentures.

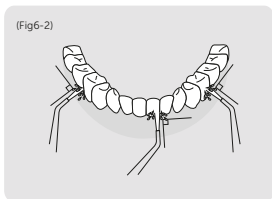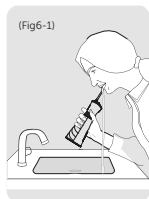Reminder!

Do not spray water directly into the periodontal pocket - it may damage the gums. (Figure 7)

After use, press the On/Off button to shut the device down. This device has a 2-minute timing function. It will automatically shut down after 2 minutes of use. If you still need to use it, please press the On/Off button. This device has a mode memory function, which automatically defaults to the last used mode when the device is turned on.

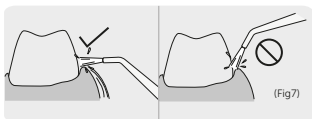
After use, do the following steps:

1. Open the water tank cover and pour out the water. (Figure 8)
2. Press the switch button to activate the water flosser and automatically drain any accumulated water. (Figure 9)
3. Press the On/Off button to shut down the water flosser, then close the water tank cover.
4. Use a cloth to dry the device. (Figure 10)

Reminder!

To avoid the spread of bacteria, ensure that there is no residual water inside the water flosser. If the device is not planned to be used for a long time, clean it and put for storage after it dries.

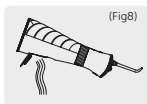Daily maintenance:

1. Clean once before each use and perform regular cleaning after use; if not used for a long time, it needs to be cleaned every six months.
2. Clean the device with pure water or neutral detergent. Do not use aggressive or abrasive detergents and hot water above 50°C for cleaning to avoid damaging the device.

Body cleaning:

- ❗ Wipe with a dry or damp cloth.
- ⊘ Do not immerse the body in water for cleaning.

Water tank cleaning:

- ❗ Wipe with a dry or damp cloth or remove the tank and rinse with pure water.
- ❗ Do not keep water inside in case the tank is not used for a long time.

Nozzle cleaning:

- ❗ Wipe the nozzle with a dry or damp cloth or remove it and rinse with pure water before wiping.
- ⊘ Do not bend, pull or twist the nozzle water pipe.

Removing the battery:

- The battery must be removed from the device before disposal.
- The device must be disconnected from the mains when removing the battery.

This symbol on the device, packaging and/or accompanying documentation indicates that used electrical and electronic products should not be disposed of with normal household waste. They should be taken to a specialized collection points

This symbol on the device, packaging, and/or accompanying documentation indicates that the battery or battery pack used in this product must not be disposed of at the end of its service life together with other household waste

Contact your local authorities for more information on existing waste collection systems. Proper disposal will help conserve valuable resources and prevent potential negative impacts on human health and the environment that may result from improper waste management

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