SCARLETT®

RUS

ENG

KAZ

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ENG

The safety measures and instructions contained in this manual do not cover all possible situations that may arise during the operation of the device. When working with the device, the user should be guided by common sense, be careful and attentive.

SECURITY MEASURES

- The manufacturer is not responsible for damages caused by non-compliance with safety regulations and rules of operation of the device.
- This device is designed for domestic use and can be used in apartments, country houses or other similar conditions of non-industrial operation. Industrial or any other inappropriate use of the device will be considered a violation of the conditions of proper operation of the device. In this case, the manufacturer is not responsible for the possible consequences.
- Do not extend the Type-C cable in doorways or near heat sources. Make sure that the cable does not twist or bend, does not come into contact with sharp objects, corners and edges of furniture.
- Improper handling may result in damage to the device, damage to property and cause harm to the health of the user.
- If the power cord is damaged, it should be replaced with a special cord obtained from the manufacturer or service center.
- The manufacturer reserves the right to make minor changes in the design of the device that do not drastically affect its safety, performance and functionality.
- Do not attempt to repair the device yourself. In case of malfunctions, please contact your dealer.
- The date of manufacture is indicated on the device and/or on the packaging, as well as in the accompanying documentation, in the format XX.XXXX, where the first two digits "XX" are the month of manufacture, the next four digits "XXXX". XX - the year of manufacture.

ATTENTION! THE WATER FLOSSER IS NOT A SUBSTITUTE FOR A TOOTHBRUSH, THEREFORE, IN ORDER TO ENSURE COMPREHENSIVE ORAL CARE OF THE ORAL CAVITY, IT IS RECOMMENDED TO MAKE BRUSHING FIRST. THE DEVICE IS DESIGNED TO CLEAN TEETH AND ORAL CAVITY, USE IT ONLY FOR ITS INTENDED PURPOSE.

REMEMBER: accidental damage to the Type-C cable can lead to problems that do not comply with the warranty conditions, as well as can cause electric shock.

- It is forbidden to operate the device outdoors: moisture or foreign objects entering the device body may cause serious damage to it. Do not charge the device near water-filled containers (bathtubs, sinks, etc.). Keep the Type-C cable away from water.

ATTENTION!

- It is not recommended to immerse the device body in water or other liquids for a long time. Before charging, make sure that the device is completely dry.
- Consult a dentist before using the device.
- When wearing dentures, implants, crowns, etc., make sure that they can be cleaned with a water flosser. The manufacturer is not responsible for possible damage.
- Do not direct a stream of water under the tongue, into the ears, nose, eyes and other sensitive areas of the body – this can lead to injury.
- Take out all removable appliances from the mouth before using the device.
- Before cleaning the device, make sure it is unplugged and turned off. Strictly follow the cleaning instructions.
- The device is not intended for use by people (including children) with reduced physical, sensory or mental abilities or with lack of experience or knowledge, unless these people are supervised or instructed in how to use this device by a person responsible for their safety.
- Children should be supervised to prevent playing with the device.
- Self-repair of the device or modification of its design is prohibited. Unprofessionally performed work can lead to damage to the device, injury and property damage.

ATTENTION! IT IS FORBIDDEN TO USE THE DEVICE IN CASE OF ANY MALFUNCTION.

SPECIFICATIONS:

Model:	SC-WF02
Device performance:	3.7 V 10W
Charging performance:	5V 1A
Battery capacity:	2000 mAh
Protection against electric shock:	Class III
The principle of operation:	pulse
Minimum water pressure:	186 kPa
Maximum water pressure:	930 kPa
Volume of the water tank:	310 ml
Jet pulsation frequency:	1400 ripples/min
Water consumption:	235 ml/min
Number of operating modes:	4
Full charge time:	up to 4 hours
Degree of moisture protection:	IPX7
Noise level:	up to 64 dB
Overall dimensions (W x H x D)	85 x 219 x 65 mm
Net weight:	275 g

PACKAGE CONTENTS

- Water Flosser
- Type-C Cable (Fig.2)
- Standard nozzles for her and for him (Fig.2)
- Orthodontic nozzle (Fig.2)
- Periodontal nozzle (Fig.2)
- Operation Manual
- Warranty card

MODEL CONFIGURATION (FIG.1)

1. Nozzle release button
2. Power button
3. Switching mode button
4. Mode indicator/power indicator
5. Water tank window
6. Replaceable nozzle
7. Installation port for the nozzle
8. Charging port
9. Water injection hole
10. Water tank cover release button
11. Water tank cover

I. BEFORE STARTING TO USE

Carefully remove the device and its components from the box.

Remove all packaging materials and promotional stickers.

Be sure to keep the warning stickers, signposting stickers (if available) and the unit serial number plate mounted to its body at their places!

Unwind the cable completely. Wipe the device body with a damp cloth. Make sure the external and visible internal parts of the device are free of.



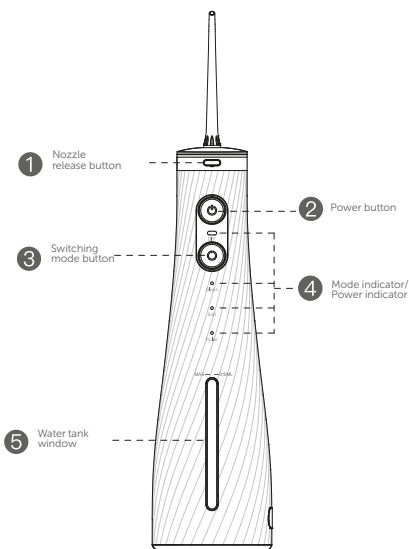
Fig.1

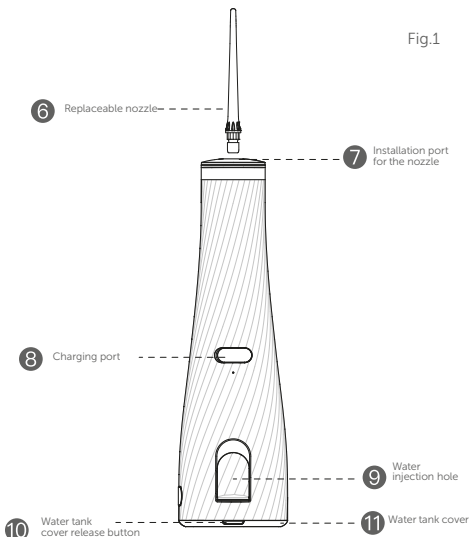ACCESSORIES

Fig.2

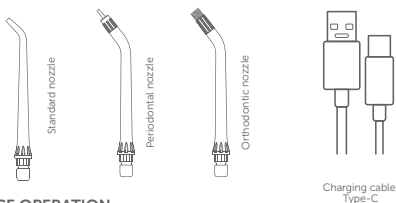II. DEVICE OPERATION

Charging the instrument

The battery should be fully charged before using the instrument. Use the Type-C cable to connect the instrument to an external power source (not included). The instrument will take up to 4 hours to fully charge.

When fully charged, the device will turn off on its own to prevent over-heating of the device. The device is ready for use.

ATTENTION! THE DEVICE WILL NOT WORK WHILE IT IS CONNECTED TO THE MAINS.

NOZZLE INSERTING AND CHANGING

The device is supplied with 4 nozzles:

- the standard nozzle (for her and for him) cleans teeth and gums from plaque and food residues, massages the gums, improving their blood circulation.
- The orthodontic nozzle is designed for cleaning braces. Replacement is recommended every 3 to 6 months;
- The periodontal nozzle is designed for cleaning and rinsing gum pockets. Replacement is recommended every 3 months. Use the device only at low power when applying the periodontal nozzle, in order to avoid damage and injury to the gingival pockets.

Make sure the device is turned off before inserting or changing the nozzles!

Insert the selected nozzle into the nozzle hole, make sure it is fixed. To remove the inserted nozzle, press the nozzle release button and pull the nozzle up (Fig.3).

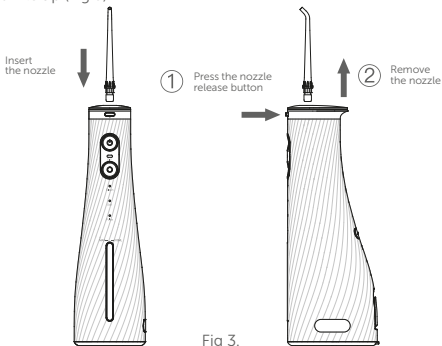
Fig 3.

USING THE DEVICE WITHOUT THE NOZZLE INSERTED IS **STRICTLY PROHIBITED!**

WORK PREPARATION:

1. Make sure the device is off.
2. Insert the selected nozzle (see the section "Nozzle inserting and changing").
3. Open the valve of the water tank and fill it no higher than the MAX mark. Close the valve (Fig 4).

To fill the tank, use purified water at room temperature or special liquids for water flossers.

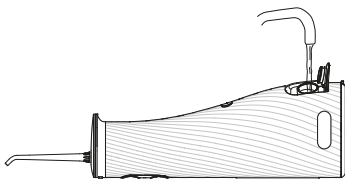
Fig.4

FILLING THE WATER TANK WITH HOT WATER, HERBAL DECOTIONS, SALT, OIL-CONTAINING OR ANTISEPTIC SOLUTIONS, TOOTHPASTE, CHEMICALS, NON-SPECIALIZED MOUTHWASH, AS WELL AS OTHER LIQUIDS NOT INTENDED FOR WATER FLOSSERS **IS STRICTLY PROHIBITED!**

4. Press the switching mode button (Fig.1) to set the desired operating mode (see the section "Operating mode switching recommendations"). The indicator of the selected mode will light up, switching is carried out in a circle. The operating modes can also be changed directly during operation of the device.
5. Press the power button (Fig.1) to activate the device.
6. Slightly open your mouth and perform the cleaning procedure directing the jet at the surface of the teeth and gums at a right angle. Move from the outermost teeth to the front ones, lingering for 2-3 seconds in the gaps between the teeth. Do not tilt the device body by more than 45° during operation. Do not block the water outlet of the nozzle! (Fig.5)
7. At the end of operation, press the power button (Fig.1), all the operation indicator will go out and the device will turn off.
8. Open the lid of the water tank (Fig.1) and drain the remaining liquid. To avoid the appearance of unpleasant odors and deposits on the walls of the tank, drain the liquid after each use of the water flosser.
9. Clean the device according to the recommendations in the "Device care" section.

OPERATING MODE SWITCHING RECOMMENDATIONS

Use the switching mode button (Fig.1)


Clean - the standard operating mode. Suitable for daily use

Soft - the soft cleaning mode. Suitable for sensitive gums

Pulse - the pulsation mode. Suitable for deep cleansing and massage of gums

DIY - the manual adjustment mode, which allows to adjust the jet pressure. Press and hold the switching mode button (Fig.1) during the DIY mode device operation - the pressure will change while the button is held.

Seven pressure levels are available, switching is carried out in a circle. When the required pressure level is reached, release the button and the device will remember the setting made. It is recommended to use the device in Soft mode during the first cleaning procedures to prevent the appearance of unpleasant or painful sensations.

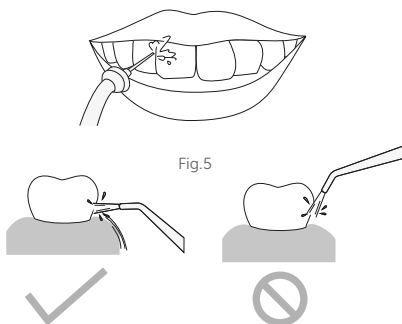III. DEVICE MAINTENANCE

Make sure the device is unplugged and turned off before cleaning. Clean the product body with a soft, damp cloth or sponge. Rinse the nozzles under running water after each use.

The inner walls of the water tank can be cleaned with a damp soft cloth if necessary (Fig. 6).

ATTENTION! IMMERSING THE DEVICE BODY IN WATER OR OTHER LIQUIDS FOR A LONG TIME IS NOT RECOMMENDED.

THE USE OF COARSE WIPES OR SPONGES, ABRASIVE PASTES WHEN CLEANING THE DEVICE **IS STRICTLY PROHIBITED.** THE USE OF ANY CHEMICALLY AGGRESSIVE SUBSTANCES IS UNACCEPTABLE.

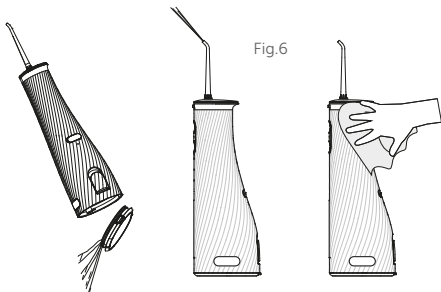
Fig.6

REMOVING THE BATTERY

- The battery must be removed from the device before disposal.
- The device must be disconnected from the mains when removing the battery.

STORAGE AND TRANSPORTATION

Store the device in a dry, ventilated place away from heating appliances and direct sunlight. It is forbidden to expose the device to mechanical shocks that may damage the device and/or violate the integrity of the package during transportation and storage. The device package should be protected from water ingress.

This symbol on the device, packaging and/or accompanying documentation indicates that used electrical and electronic products should not be disposed of with normal household waste. They should be taken to a specialized collection points

This symbol on the device, packaging, and/or accompanying documentation indicates that the battery or battery pack used in this product must not be disposed of at the end of its service life together with other household waste

Contact your local authorities for more information on existing waste collection systems. Proper disposal will help conserve valuable resources and prevent potential negative impacts on human health and the environment that may result from improper waste management

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